allard

Instructions For Use

Kiddie Composite Footplate

Art. no. 28037

Instructions for Use / Användarinstruktion
Brugsanvisning / Anvendelseinstruktion
Käyttöohjeet / Gebrauchsanweisung
Gebruiksaanwijzing / Instructions d'utilisation
Istruzioni per l'uso / Instrucciones de uso
Instruções de Uso / Kasutusjuhend
Lietošanas instrukcija / Naudojimo instrukcija

Product description

A composite footplate that limits the range-of motion (ROM) in the meta-tarso-phalangeal joints (mtp-joints), distributes the pressure in the forefoot and facilitates the roll-over.

Intended use

The composite foot plate is intended to protect the meta-tarso-phalangeal joints (mtp-joints) and the forefoot by limiting the range-of motion (ROM), distributing the pressure in the forefoot and facilitating roll-over. The composite foot plate can also be used after immobilizing treatment with POST-OP shoe or walker, in order to gradually increase movement. It can also be used to improve balance after toe amputation, dig IV and/or V. Must be covered by an insole. Single patient Multiple use.

Indications

Toe fractures, stress fractures, and other injuries or arthritis conditions in the forefoot/midfoot, forefoot trauma, toe/mid-foot sprain, post op recovery, toe amputation dig IV-V, forefoot dysmelia, Idiopathic toe-walkers

Contraindications

Open ulcers on the foot. Non-stable fractures.

Adjustments

Fitting and modification of this product may only be performed by qualified healthcare professionals.

Warnings

Extra precaution needed if the patient presents with conditions that may lead to ulcers . Make sure that the composite footplate is covered by a proper and functional foot orthotic that manages the pressure distribution on the plantar surface and that the selected shoe has sufficient inner volume to ensure that pressure sores do not occur.

Fiber reinforced thermoset composite: Do not heat the material. When grinding/cutting: Protect the eyes, use dust protection over nose and mouth, make sure that the edges of the footplate are carefully smoothed and all loose fiber particles are removed. Do not overheat (max. 100° C/212° F) when grinding the composite.

Material Specifications

Carbon fibre in an epoxy matrix

User instruction

Check the condition of your foot daily. In the case of skin irritation or discoloration, discontinue use immediately and contact your fitter/prescriber. If your ability to monitor the foot is reduced, ask for assistance. Do not bend the composite footplate excessively, for example, by squatting.

Washing Instructions

Wash with soap and water.

Questions

For any questions or issues, please contact your healthcare professional/fitter.

Disposal

Disposal in accordance with local legislation.

Information

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Fitting the composite footplate

1


Remove the insole from the shoe. Choose the composite footplate closest to the patient's shoe size.

2

Use the insole from the footwear as a template. Trace the template onto the composite footplate. The insole should be a couple of millimeters larger than the composite footplate. Grind the composite footplate with a belt sander or Dremel and smooth the edges.

3

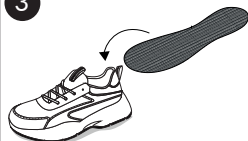Place the composite footplate in the shoe.

4

Return the shoe insole on top of the composite footplate.

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©Camp Scandinavia_0626_Rev. 01
Camp Scandinavia AB, Karbingatan 38