

MedSpec Product Range



Explore the Wrist & Elbow Range

Stability and comfort for hands and elbows



Wrist Lacer™ II, Short

Highly stable wrist orthosis with mouldable aluminium splints on both the dorsal and volar sides. The splints can be removed and shaped as needed. A cushion on the inside of the orthosis sleeve protects the tendons at the base of the thumb near the distal lacing.

Made from durable suede-like material with a breathable lining. Locking system with a unique design that provides even compression across all four straps. Compression can be easily adjusted for patients with atypical forearm and wrist shapes. One-handed operation system. Adjustable thumb grip.

Recommended Range Of Application: Carpal tunnel syndrome, post-traumatic, and post-operative use.



Wrist Lacer™ II, Long

Highly stable long wrist orthosis with removable and mouldable aluminium splints on both dorsal and volar sides. A cushion on the inside of the orthosis sleeve protects the tendons at the base of the thumb near the distal lacing.

Made from durable suede-like material with a breathable lining. Unique locking system providing even compression across all four straps. Compression can be easily adjusted for patients with atypical forearm or wrist shapes. One-handed lacing system. Adjustable thumb support.

Recommended Range Of Application: Post-traumatic, post-operative, rheumatologic conditions and spasticity.



Tripod™ Wrist Lacer II.

The Tripod Wrist Lacer II is a supportive wrist brace designed for enhanced immobilization and easy application. Ulnar, palmar, and dorsal aluminum stays minimize wrist flexion-extension and also help control radial-ulnar deviation for effective stabilization.

A patented lace closure design distributes equal tension across all laces, creating uniform compression for outstanding support. The tension can be easily adjusted for patients with

disproportionate forearm and wrist circumferences. Made from a durable polymer felt liner and perforated suede that allows moisture to evaporate, enhancing comfort and patient compliance.

Recommended Range Of Application: Post-fracture care, post-operative care, carpal tunnel syndrome, wrist fractures, or wrist strain.



Ryno Lacer® II wrist orthosis

Highly stable wrist and thumb orthosis with mouldable aluminium splints on the dorsal, volar, and radial sides. The radial thumb splint extends along the full length of the splint and effectively immobilises the thumb.

Made from durable suede-like material with a breathable lining. The lacing locking system is easy to open and close with one hand.

Recommended Range Of Application: Post-traumatic, post-operative, rheumatologic conditions and spasticity.



Cubital Tunnel Brace

The Cubital Tunnel Brace is a lowprofile night splint that comfortably holds the elbow in a neutral position for the treatment of Cubital Tunnel Syndrome.

Ventral design, with no exposed rigid materials. Enhances patient compliance during sleep. This Cubital Tunnel Brace features an adjustable aluminum stay that allows the elbow to be positioned at the optimum angle.

Unique strap retention system allows for easy application by holding the proximal and distal closure straps within their plastic buckles. The middle strap may be removed if desired for various treatment protocols.

Recommended Range Of Application: Cubital tunnel syndrome. Ulnar nerve entrapment at the elbow.

Explore the Foot & Ankle Range

Support, stability and comfort for every step



Supinator™ PTT Stabiliser

The Supinator is a comfortable and relatively low-profile brace for the treatment of posterior Tibial Tendon Dysfunction (PTTD). With its patent-pending design, the Supinator reduces over-pronation by applying a supinating force.

The arch strap works together with a repositionable arch pad to lift the medial arch and the heel strap restricts eversion of the calcaneus (heel). The arch pad may be removed if desired to allow the Supinator PTTD brace to be worn with an orthotic. In addition, the Supinator features the Speed Lacer II closure system for ease of application and equal tension across all laces.

Recommended Range Of Application:
Posterior tibial tendon dysfunction (PTTD), acquired adult flatfoot deformity posterior tibial tendon insufficiency.

Contraindications:
Lateral ankle sprain, lateral ankle instability.



ASO® Ankle Stabiliser

Stable low-profile ankle support with front lacing. The rigid nylon straps are positioned to mimic taping and prevent the foot from supinating or pronating.

The straps can be adjusted depending on the nature of the injury. For added comfort, the orthosis features new reinforced Coolflex™ padding along the inside, next to the Achilles Tendon.

Recommended Range Of Application:
Acute or chronic ankle instability, and for preventive use during sports activities.



EVO™ Ankle Stabiliser

Stable low-profile ankle support with stabilising front lacing. The rigid nylon straps mimic taping and prevent the foot from supinating or pronating. On the ASO EVO, the nylon straps around the calcaneus provide greater resistance to calcaneal inversion/ eversion than the standard ASO. The upper part is reinforced with a plastic cuff, enhancing protection against high ankle sprains.

For added comfort, the orthosis features new reinforced Coolflex™ padding along the inside, next to the Achilles tendon.


Recommended Range Of Application:
Acute or chronic ankle instability, preventive use during sports activities, high ankle sprains, and subtalar joint instability.

ASO® & EVO™ Sizing Charts

Use your shoe size as a guide only. For the best fit, measure your ankle circumferences as shown below .

Ankle Circumference

Paediatric	23–25 cm
X-Small	25–28 cm
Small	28–30 cm
Medium	30–33 cm
Large	33–38 cm
X-Large	38–41 cm
XX-Large	41–43 cm
XXX-Large	43+ cm

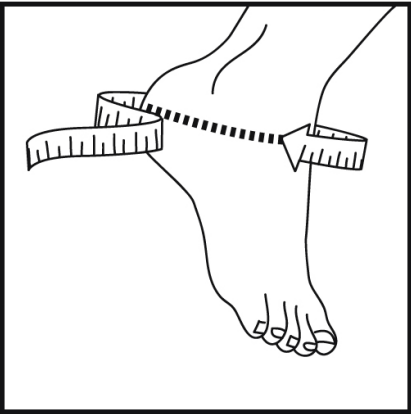

Men's shoe size chart

Width	15–18 cm	19–23 cm	24–28 cm	29–33 cm	34–38 cm	38 cm
Narrow	Paediatric	X-Small	Small	Medium	Large	X-Large
Regular	X-Small	Small	Medium	Large	X-Large	XX-Large
Wide	Small	Medium	Large	X-Large	XX-Large	XXX-Large

Women's shoe size chart

Width	13–15 cm	17–20 cm	22–25 cm	27–30 cm	32–36 cm
Narrow	Paediatric	X-Small	Small	Medium	Large
Regular	X-Small	Small	Medium	Large	X-Large
Wide	Small	Medium	Large	X-Large	XX-Large

Ordering Information

Wrist Lacer™ II, Short

Item.no	Model	Size	Measurement: Wrist Circumference	L/R	Colour
285781009	Wrist Lacer II, Short	XX-Small	<12 cm	Left	Black
285781010	Wrist Lacer II, Short	X-Small	12-14 cm	Left	Black
285781011	Wrist Lacer II, Short	Small	14-16 cm	Left	Black
285781012	Wrist Lacer II, Short	Medium	16-18 cm	Left	Black
285781013	Wrist Lacer II, Short	Large	18-20 cm	Left	Black
285781014	Wrist Lacer II, Short	X-Large	20-22 cm	Left	Black
285781015	Wrist Lacer II, Short	XX-Large	22-25 cm	Left	Black
285782009	Wrist Lacer II, Short	XX-Small	<12 cm	Right	Black
285782010	Wrist Lacer II, Short	X-Small	12-14 cm	Right	Black
285782011	Wrist Lacer II, Short	Small	14-16 cm	Right	Black
285782012	Wrist Lacer II, Short	Medium	16-18 cm	Right	Black
285782013	Wrist Lacer II, Short	Large	18-20 cm	Right	Black
285782014	Wrist Lacer II, Short	X-Large	20-22 cm	Right	Black
285782015	Wrist Lacer II, Short	XX-Large	22-25 cm	Right	Black

Wrist Lacer™ II, Long

Item.no	Model	Size	Measurement: Wrist Circumference	L/R	Colour
285791010	Wrist Lacer II, Long	X-Small	12-14 cm	Left	Black
285791011	Wrist Lacer II, Long	Small	14-16 cm	Left	Black
285791012	Wrist Lacer II, Long	Medium	16-18 cm	Left	Black
285791013	Wrist Lacer II, Long	Large	18-20 cm	Left	Black
285791014	Wrist Lacer II, Long	X-Large	20-22 cm	Left	Black
285792010	Wrist Lacer II, Long	X-Small	12-14 cm	Right	Black
285792011	Wrist Lacer II, Long	Small	14-16 cm	Right	Black
285792012	Wrist Lacer II, Long	Medium	16-18 cm	Right	Black
285792013	Wrist Lacer II, Long	Large	18-20 cm	Right	Black
285792014	Wrist Lacer II, Long	X-Large	20-22 cm	Right	Black

Tripod™ Wrist Lacer II

Item.no	Model	Size	Measurement: Wrist Circumference	L/R	Colour
285811010	Tripod™ Wrist Lacer II	X-Small	12-14 cm	Left	Black
285811011	Tripod™ Wrist Lacer II	Small	14-16 cm	Left	Black
285811012	Tripod™ Wrist Lacer II	Medium	16-18 cm	Left	Black
285811013	Tripod™ Wrist Lacer II	Large	18-20 cm	Left	Black
285811014	Tripod™ Wrist Lacer II	X-Large	20-22 cm	Left	Black
285812010	Tripod™ Wrist Lacer II	X-Small	12-14 cm	Right	Black
285812011	Tripod™ Wrist Lacer II	Small	14-16 cm	Right	Black
285812012	Tripod™ Wrist Lacer II	Medium	16-18 cm	Right	Black
285812013	Tripod™ Wrist Lacer II	Large	18-20 cm	Right	Black
285812014	Tripod™ Wrist Lacer II	X-Large	20-22 cm	Right	Black

Supinator™ PTT Stabiliser

Item.no	Model	Size	Measurement: Ankle Circumference	L/R	Colour
281961010	Supinator™ PTT Stabiliser	X-Small	25,4–27,9 cm	Left	Black
281961011	Supinator™ PTT Stabiliser	Small	27,9–30,5 cm	Left	Black
281961012	Supinator™ PTT Stabiliser	Medium	30,5–33 cm	Left	Black
281961013	Supinator™ PTT Stabiliser	Large	33–35,6 cm	Left	Black
281961014	Supinator™ PTT Stabiliser	X-Large	35,6–38,1 cm	Left	Black
281962010	Supinator™ PTT Stabiliser	X-Small	25,4–27,9 cm	Right	Black
281962011	Supinator™ PTT Stabiliser	Small	27,9–30,5 cm	Right	Black
281962012	Supinator™ PTT Stabiliser	Medium	30,5–33 cm	Right	Black
281962013	Supinator™ PTT Stabiliser	Large	33–35,6 cm	Right	Black
281962014	Supinator™ PTT Stabiliser	X-Large	35,6–38,1 cm	Right	Black

ASO® Ankle Stabiliser

Item.no	Model	Size	Measurement: Ankle Circumference	L/R	Colour
281860010	ASO® Ankle Stabiliser w/Stays	X-Small	25-27,5 cm	Bilateral	Black
281860011	ASO® Ankle Stabiliser w/Stays	Small	27,5-30 cm	Bilateral	Black
281860012	ASO® Ankle Stabiliser w/Stays	Medium	30-32,5 cm	Bilateral	Black
281860013	ASO® Ankle Stabiliser w/Stays	Large	32,5-35 cm	Bilateral	Black
281860014	ASO® Ankle Stabiliser w/Stays	X-Large	35-37,5 cm	Bilateral	Black
281860015	ASO® Ankle Stabiliser w/Stays	XX-Large	37,5-40 cm	Bilateral	Black
281880009	ASO® Ankle Stabiliser	XX-Small	22,5-25 cm	Bilateral	Black
281880010	ASO® Ankle Stabiliser	X-Small	25-27,5 cm	Bilateral	Black
281880011	ASO® Ankle Stabiliser	Small	27,5-30 cm	Bilateral	Black
281880012	ASO® Ankle Stabiliser	Medium	30-32,5 cm	Bilateral	Black
281880013	ASO® Ankle Stabiliser	Large	32,5-35 cm	Bilateral	Black
281880014	ASO® Ankle Stabiliser	X-Large	35-37,5 cm	Bilateral	Black
281880015	ASO® Ankle Stabiliser	XX-Large	37,5-40 cm	Bilateral	Black
281880016	ASO® Ankle Stabiliser	XXX-Large	>40 cm	Bilateral	Black

EVO™ Ankle Stabiliser

Item.no	Model	Size	Measurement: Ankle Circumference	L/R	Colour
281890010	EVO™ Ankle Stabiliser	X-Small	25-27,5 cm	Bilateral	Black
281890011	EVO™ Ankle Stabiliser	Small	27,5-30 cm	Bilateral	Black
281890012	EVO™ Ankle Stabiliser	Medium	30-32,5 cm	Bilateral	Black
281890013	EVO™ Ankle Stabiliser	Large	32,5-35 cm	Bilateral	Black
281890014	EVO™ Ankle Stabiliser	X-Large	35-37,5 cm	Bilateral	Black
281890015	EVO™ Ankle Stabiliser	XX-Large	37,5-40 cm	Bilateral	Black



Support for Better Life

Everyone should be able to live their life to the fullest, regardless of their mobility challenges. With innovative solutions developed in close collaboration with healthcare professionals and patients, we strive to provide Support for Better Life.

allard

Allard Support UK Ltd.

Coes Road, Dundalk, Co. Louth, A91 E868 | +353 42 932 8177 | info@allarduk.co.uk | allarduk.co.uk