

Selection[®] Ulnar Deviation Splint

- Low Profile Design, meets the shape and contours of the hand
- Breathable fabric
- Fabric can be trimmed
- Easy fastening
- Can be worn with Selection[®] Wrist Splint
- Machine washable



Improve finger alignment to enhance hand function and functional grip

Selection® Ulnar Deviation Splint is intended to improve finger alignment and enhance hand function and functional grip.

This hand splint is exclusively designed and produced using our unique three layer Selection fabric. With its low profile design, it meets the contours of the hand without impeding function. The fabric can be trimmed without fraying, and therefore allows further customization if required. The easy fastening surface allows the hook to be placed anywhere on the splint to optimise the directional pull needed.

The Ulnar Deviation splint corrects alignment of joints of the fingers to improve extension, strength and function of the hand and thereby can provide reduced pain and joint stiffness. The orthosis is designed so that the material extends far to the ulnar side at the 5th MCP, to support the best possible direction of pull and distribute the pressure.

Should you need to provide wrist stabilisation for your patient, this splint can also be worn in combination with a Selection® wrist splint.

Intended use

Ulnar deviation is likely to be a characteristic for those with various forms of inflammatory disease, such as Rheumatoid Arthritis. The Ulnar Deviation Splint may also be used for realignment and stabilization in patients with joint laxity and muscle weakness. Not intended for use with fixed contractures in the hand.



Item No. Left	Item No. Right	Size	Circumference MCP 2-5 (cm)
35218 1010	35218 2010	X-Small	13-16 cm
35218 1011	35218 2011	Small	16-19 cm
35218 1012	35218 2012	Medium	19-22 cm
35218 1013	35218 2013	Large	22-25 cm

Developed for optimal function and comfort

Selection® - a full range of hand and wrist orthoses

Our Selection range is created for anatomically compliant orthoses that provide support and stability where needed, but at the same time are flexible and allow movement so that the body part can be functional. Our mission with Selection is for the products to feel like a second skin and to be customisable for a perfect fit for everyone.


The soft interior and edges, together with the breathable and quick drying material, provides a comfortable and cool fit. Selection orthoses has lacing for customized size and support and can be cut for further customization without fraying. The low profile of the products allows them to be used under clothes and gloves.

Tested to be used close to the body according to SGS.



The Selection Material

All our materials are carefully tested at internationally approved testing institutes to ensure that they contain no substances that can be harmful to either adults or children. The Velcro straps used in the products are carefully selected to suit their particular function and has also been tested to ensure long lasting durability.

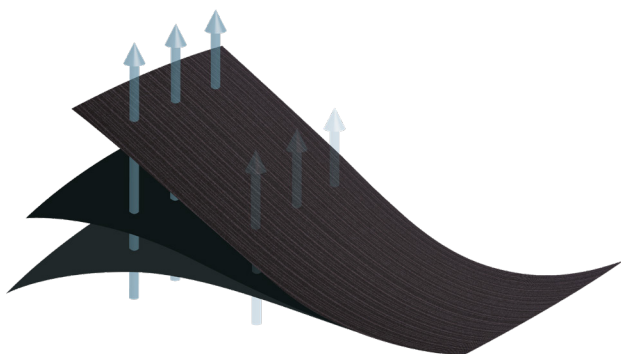

Selection Powernet

The function carrier, the inner textile mesh structure, gives the material its perfect elasticity and recovery for good consistent support.

Selection Air

Thin, soft and 100% air-permeable material. Used in selected parts of the orthosis where the textile should feel invisible and not retain heat. Gives extra smooth function.

An excellent combination

The combination provides an anatomical orthosis that allows unrestricted movement and stability. A perfect fit that is soft inside and durable outside. Elegant with high comfort.

Support for Better Life

Everyone should be able to live their life to the fullest, regardless of their mobility challenges. With innovative solutions developed in close collaboration with healthcare professionals and patients, we strive to provide Support for Better Life.

allard

Allard Support UK Ltd.

Coes Road, Dundalk, Co. Louth, A91 E868 | +353 42 932 8177 | info@allarduk.co.uk | allarduk.co.uk