

allard|_{UK}

SO AIRY
& LIGHT
you barely feel it!



X-LITE® CLASSIC - NOW IN NEW COLORS



We are happy to present our well-known X-LITE® CLASSIC material in two new attractive colors. Anthracite and Royal Blue.

- X-LITE® is an airy and lightweight fabrication material engineered to meet your splinting, bracing and casting demands.
- X-LITE® is made from a 100% cotton mesh that has been impregnated with a non-toxic thermoplastic resin. All reactive substances in the thermoplastic have been cured, which means it does not contain any residues in unreacted form.
- X-LITE® is a NON-Toxic material and emits no harmful substances when heated, applied or in the drying process. The material can be applied without the use of gloves and special ventilation equipment in the casting room. It is safe for all users to work with X-LITE®, including pregnant woman, asthmatic patients and people that suffer from allergic reactions.
- X-LITE® meets the international standard of biodegradation and is sorted as non-hazardous waste. A true environmentally friendly material!
- X-LITE® is heated in a waterbath set at 70°C.
- Material starts to set after 1 min. Fully cured after 20 min.

X-LITE® is a NON-Toxic and environmentally friendly product providing a safe working environment and protecting our nature.




100%
Natural Cotton

+




Non-Toxic
Polymer

=

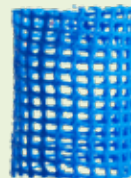

X-LITE®

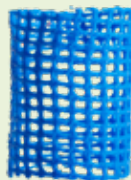

X-LITE®

+



Compost

=



CO₂

&



Water



Our environment - our future - your choice!

FEATURES AND BENEFITS X-LITE® CLASSIC

- Ultra-light providing an airy and breathable splint
- Strong yet low profile
- Moldable after only 10 seconds in the Waterbath
- Add layers for increased strength
- Reheat any number of times
- Remove water and mould directly on the patient
- No fingerprints
- Hook, padding and other devices are easily attached for a permanent bond



Item No.	Width	Length	Colour	Package
64136 0000	7,5 cm	1,8 m	Anthracite	5 rolls
64146 0000	10 cm	1,8 m	Anthracite	5 rolls
64168 0000	15 cm	1,8 m	Anthracite	2 rolls
64302 0000	20 cm	37,5 cm	Anthracite	8 sheets
64303 0000	45 cm	50 cm	Anthracite	4 sheets
64255 0000	15 cm	10 m	Anthracite	1 dispenser
64137 0000	7,5 cm	1,8 m	Royal Blue	5 rolls
64147 0000	10 cm	1,8 m	Royal Blue	5 rolls
64169 0000	15 cm	1,8 m	Royal Blue	2 rolls
64304 0000	20 cm	37,5 cm	Royal Blue	8 sheets
64305 0000	45 cm	50 cm	Royal Blue	4 sheets
64256 0000	15 cm	10 m	Royal Blue	1 dispenser


Our environment - our future - your choice!

SO AIRY & LIGHT you barely feel it!

X-LITE® PLUS EDGING TAPE

The Edging Tape gives a neat finish to the splint and secures a soft edge. With a back of Premium material it provides extra stability to the splint. Works well with all our materials including 3D-LITE™. Heat in 70°C waterbath.

Item No.	Colour	Width	Length	Package
64235 0000	White	25 mm	4 m	Dispenser box
64236 0000	Red	25 mm	4 m	Dispenser box
64237 0000	Blue	25 mm	4 m	Dispenser box
64238 0000	Black	25 mm	4 m	Dispenser box
64239 0000	Hot Pink	25 mm	4 m	Dispenser box

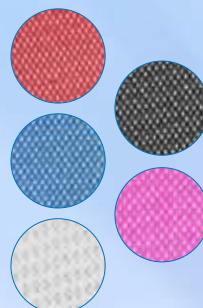

X-LITE® THERMOPLASTIC HOOK

An easy solution to attach your pile closure straps to your splint or brace. X-LITE® polymer is on the back of the hook and bonds instantly and permanently on all our materials. Heat the back with the heatgun or in the waterbath – spot heat the splint for a permanent bond.

Item No.	Colour	Width	Length	Package
64195 0000	White	50 mm	50 cm	5 Sheets
64196 0000	White	50 mm	5 m	1 Dispenser
64197 0000	Black	50 mm	50 cm	5 Sheets
64198 0000	Black	50 mm	5 m	1 Dispenser

NEW!
AVAILABLE IN BLACK!



allard | UK

www.allarduk.co.uk

Allard Support UK, Ltd.
Coes Road Dundalk A91 E868
Co. Louth, Ireland

Phone: +353 42 932 8177
Fax: +353 42 932 8182
FreePhone: 0800 051 7061
FreeFax: 0800 051 7069
customerservice@allarduk.co.uk